

SINCE 1925

**KEYSTONE**

SAFETY ↔ CLEANROOM PRODUCTS

## BLUE BHD STERILIZATION WRAPPER



When secured with autoclave tape, or elastic bands, Blue BHD wrappers provide an effective microbial barrier from external sources of contamination.

- Treated to impart water and alcohol repellency
- Superior bacterial barrier properties compared to linen
- Sterilization indicator tapes adhere firmly and remove easily
- Ideal choice for for light trays, small basins and towels
- Resists tearing and punctures
- Does not contain any natural latex rubber
- Processed and packaged in ISO certified cleanrooms
- Functional equivalent to the discontinued Cardinal Health™ Bio-Shield® wrap

These versatile wraps offer fluid repellency and provide an excellent barrier against airborne bacteria. The soft, drapable non-woven blend of cellulose and polyester is virtually lint-free and helps protect wrapped contents from contamination.

These wraps are wet formed, lightweight non-woven fabric comprised of natural wood pulp and polyester fibers bonded with a resin binder and are designed for steam sterilization; but also compatible with EtO or gamma irradiation.

### MATERIAL SPECIFICATIONS:


PROPERTY	UNITS	BHD MATERIAL
Basis Weight	Gm/m2	57
Dry Tensile Strength (MD)	grams	8,100
Dry Tensile Strength (CD)	grams	7,100
Elongation Dry	%	27.3
Air Permeability	liters / 100 sq. cm	110
Suter Water Head (pressure test)	mm	400
Dry Mullen Burst	Grams / sq.cm	2,150
Flammability	seconds	8

SINCE 1925

**KEYSTONE**

SAFETY ↔ CLEANROOM PRODUCTS

# BLUE BHD STERILIZATION WRAPPER

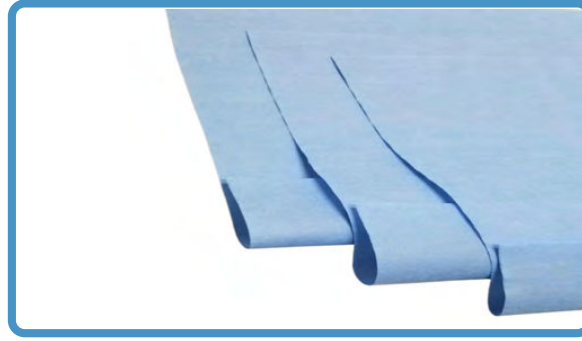

KEYSTONE PART NUMBER	DIMENSIONS	KEYSTONE PART NUMBER	DIMENSIONS
BHDSHEET0303	3" X 3"	BHDSHEET2020	20" X 20"
BHDSHEET0404	4" X 4"	BHDSHEET2424	24" X 24"
BHDSHEET0606	6" X 6"	BHDSHEET3030	30" X 30"
BHDSHEET0808	8" X 8"	BHDSHEET3636	36" X 36"
BHDSHEET0909	9" X 9"	BHDSHEET4040	40" X 40"
BHDSHEET0911	9" X 11"	BHDSHEET4545	45" X 45"
BHDSHEET1212	12" X 12"	BHDSHEET5454	54" X 54"
BHDSHEET1515	15" X 15"	BHDSHEET5472	54" X 72"
BHDSHEET1818	18" X 18"	BHDSHEET6060	60" X 60"